

ORTHOMITE LC

LIGHT-CURE ORTHODONTIC ADHESIVE

Easy to Use, Time-saving
Reliable Bonding Performance
Applicable to Any Brackets




Orthomite LC is a light-cure adhesive resin for orthodontic Bonding of metal, ceramic and plastic brackets to tooth surface.

Easy to Use, Time-Saving

Through the innovative design of the base monomer, Orthomite LC paste can bond to etched enamel well without priming.

Reliable Bonding Performance

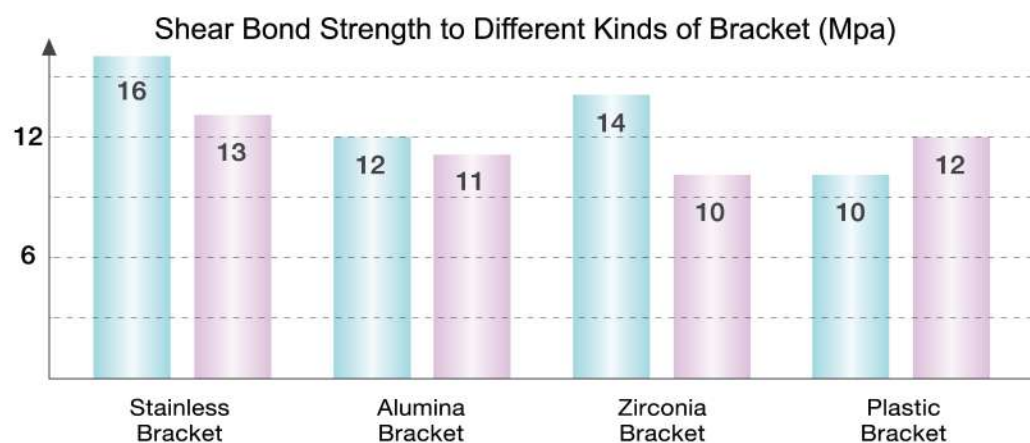

- Immersed in 370C water for 24 hrs
- Thermocycled 5,000 times (5-550C)

Applicable to any brackets
Minimal bond failure

Clinical Procedure



Remove excess adhesive Before curing.

Metal brackets: both interproximal sides 10 sec. Ceramic brackets and plastic brackets: 10 sec. through brackets.

Contents of the Kit

- | | | |
|-------------------------|---|--------|
| 1. Orthomite LC Paste | - | 4g x 2 |
| 2. Orthomite LC Etchant | - | 5mL |
| 3. Sponge (L.S) | - | 1 |